

10X Nikon CFI60 TU Plan Epi Infinity Corrected Obj

See More by [Nikon](#)



#58-516: 10XNikon CFI60 TU Plan Epi

Stock **#58-516** **8 In Stock**

-

1

+

£1,128⁰⁰

ADD TO CART

Volume Pricing	
Qty 1+	£1,128.00 each
Need More?	Request Quote

 Prices shown are exclusive of VAT/local taxes

Product Downloads

SPECIFICATIONS

General

Model Number:

MJE12100

Compatible Tube Lens Focal Length (mm):

Focal Length: 200mm

Microscope Objective	Type:
Infinity Corrected	Style:
Nikon	Manufacturer:
Physical & Mechanical Properties	
2.20	Field of View (mm):
42.40	Length excluding Threads (mm):
29	Maximum Diameter (mm):
125.00	Weight (g):
Optical Properties	
0.64mm	Horizontal Field of View, 1/2" Sensor:
0.88mm	Horizontal Field of View, 2/3" Sensor:
20.00	Focal Length FL (mm):
10X	Magnification:
0.30	Numerical Aperture NA:
1.0	Resolving Power (μm):
3.06	Depth of Field (μm):
2.5	Field of View, 25 Diameter Field Eyepiece (mm):
2.2	Field of View, 22 Diameter Field Eyepiece (mm):
17.5	Working Distance (mm):
22	Field Number (mm):
60	Parfocal Length (mm):
N/A	Immersion Liquid:
Sensor	
2/3"	Maximum Sensor Format:
Threading & Mounting	
M25 x0.75	Mounting Threads:
Regulatory Compliance	
Compliant	RoHS 2015:
View	Certificate of Conformance:

PRODUCT DETAILS



- Excellent Color Reproduction
- Long Working Distance and High NA
- High Contrast with Minimal Flare
- Strain Free

Nikon CFI60 Infinity Corrected Brightfield Objectives are designed to meet advanced imaging requirements and are the culmination of Nikon's optical technology. The objectives offer increased brightness and the highest possible combination of working distance and numerical aperture (for light gathering ability) due to the 60mm parfocal distance and a large barrel diameter. Designed to correct lateral and axial chromatic aberrations over the entire field of view, these objectives produce high contrast and high-resolution images. Nikon CFI60 Infinity Corrected Brightfield Objectives have excellent color reproduction, strain free. These brightfield objectives are offered in magnifications from 5X up to 100X.

The accessory tube lens ([#58-520](#)) can be used to focus light from the infinity corrected objective onto the image plane. For optimal performance, the distance between the tube lens and the objective's shoulder should be kept within 100-200mm. This spacing allows accessories to be inserted in the optical path to accommodate specific needs. This tube lens can be integrated into C-Mbunt systems using the Nikon Tube Lens to C-Mount Holder ([#11-150](#)).

Note: The four socket head cap screws must be removed to place the tube lens within the Nikon Tube Lens to C-Mount Holder.

TECHNICAL INFORMATION

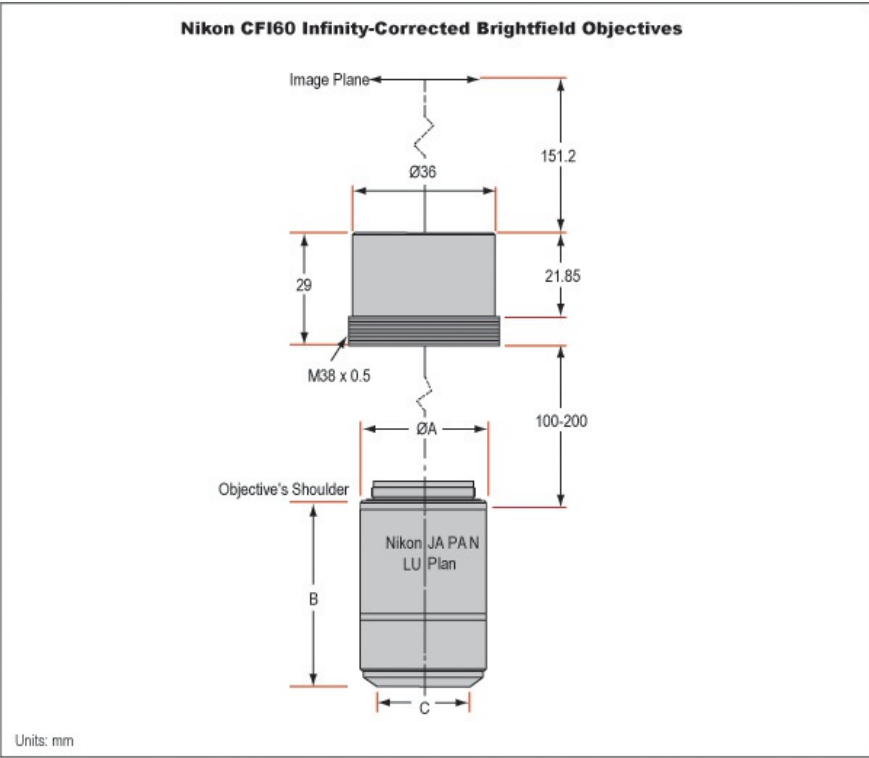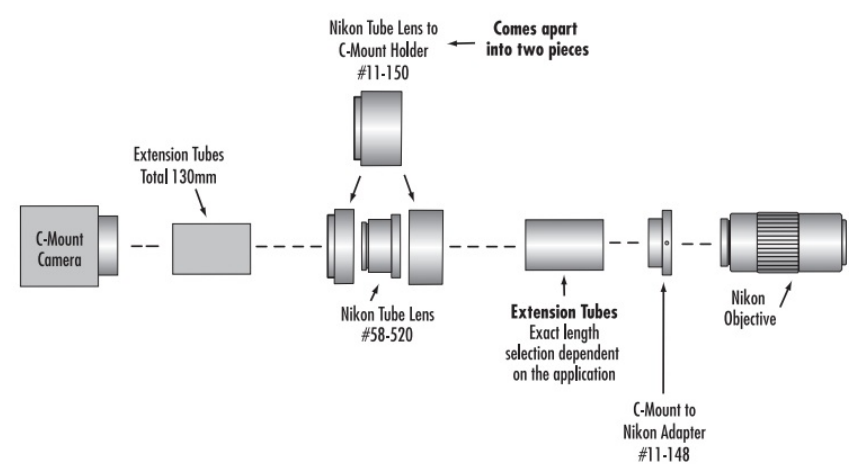

	5X	10X	20X	50X	100X
A	29.0mm	29.0mm	29.0mm	30.0mm	30.0mm
B	36.3mm	42.4mm	40.9mm	48.9mm	55.4mm
C	20.0mm	17.0mm	22.5mm	21.5mm	21.5mm

Nikon CF160 Infinity Corrected Brightfield Objectives

