

[See all 20 Products in Family](#)

## TECHSPEC® Fluorescence Filter Set for Cy3 Fluorescence Dye



Fluorescence Filter Sets



Stock #67-006 [CONTACT US](#)

1 **£760<sup>00</sup>**

**ADD TO CART**

### Volume Pricing

Qty 1-9	<b>£760.00</b> each
Qty 10-25	<b>£684.00</b> each
Qty 26-49	<b>£616.00</b> each
Need More?	<a href="#">Request Quote</a>

**!** Prices shown are exclusive of VAT/local taxes

### Product Downloads

### General

Fluorescence Filter Kit **Type:**  
CY3 **Compatible Fluorophore:**

#67-031      Excitation Filter:

#67-034      Emission Filter:

#67-082      Dichroic Filter:

## Optical Properties

**Substrate:**  Fused Silica (Corning 7980)

**Coating:** Hard Coated

**Excitation Wavelength (nm):** 513 - 556

**Emission Wavelength (nm):** 570 - 613

**Dichroic Cut-On Wavelength (nm):** 562.00

**Wavelength Range (nm):** 513 - 613

## Regulatory Compliance

**Certificate of Conformance:** [View](#)

## Product Details

The TECHSPEC® Fluorescence Filter Set for Cy3 Fluorescence Dye features an excitation filter, an emission filter, and a dichroic filter that complement the Cy3 spectrum. The excitation filter passes wavelengths within the Cy3 excitation/absorption spectrum with a range of 513 - 556nm while the emission filter passes wavelengths within the Cy3 emission spectrum with a range of 570 - 613nm.


- Excitation, Emission, and Dichroic Filters for Fluorescence Imaging
- >93% Transmission and OD 6 Blocking for Maximum Brightness and Contrast
- Hard Sputtered Coatings on Single Substrate
- [TECHSPEC Pre-Mounted Filter Cubes](#) Also Available

Our TECHSPEC® Fluorescence Filter Sets are designed to provide brilliant images when integrated into a fluorescence microscope. The high transmission of the bandpass filters guarantees the brightest possible images, while the deep blocking ensures no unwanted light gets to the detector. This ensures images will be brighter and will have much greater signal-to-noise than what is commonly available with soft coated, evaporated filters. TECHSPEC Fluorescence Filter Sets are ideal for fluorescence microscopy or fluorescence imaging applications and are available for common fluorescent dyes and fluorescent proteins. Additionally these filters are also used in polymerase chain reaction (PCR) diagnostic instruments.

Fluorescence filter sets include an Excitation and Emission filter from our [TECHSPEC® Fluorescence Bandpass Filters](#) and a Dichroic filter from our [TECHSPEC® Dichroic Filters](#). Excitation filters are placed within the illumination path of fluorescence microscopes to pass only wavelengths within the excitation range of the fluorophore under inspection. Emission filters are placed within the imaging path of a fluorescence microscope to pass only wavelengths within the emission range of the fluorophore. Dichroic filters are placed in between the excitation filter and the emission filter to reflect the excitation signal towards the fluorophore while transmitting the emission signal to the detector.

Filter sets are listed with the most common fluorophores. For a full list of compatible fluorophores with each set, view [Fluorophores and Optical Filters for Fluorescence Microscopy](#). TECHSPEC Fluorescence Filter Sets are available pre-mounted in [Olympus or Nikon Filter Cubes](#).

## Technical Information

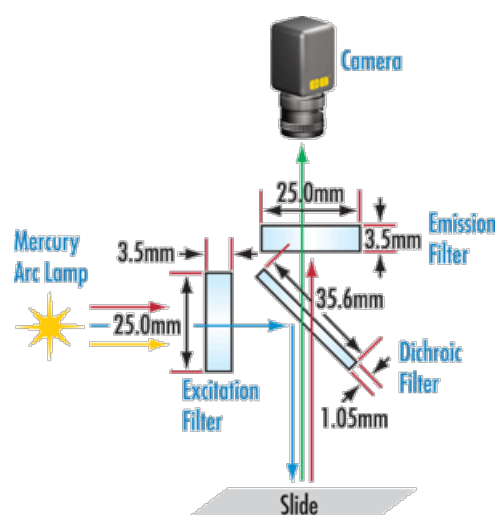

	Excitation and Emission Filters	Dichroic Filter
Size	25.0mm	25.2 x 35.6mm
Thickness	3.5mm	1.05mm
Tolerances		
Size	+0.0/-0.1	
Thickness	±0.1mm	
Clear Aperture	22.0mm	22.0 x 30.0mm
Transmitted Wavefront (RMS)	<λ/4	
Angle of Incidence	0.0° ±5.0°	45.0° ±5.0°
Compatible Microscopes	Most common Nikon, Olympus, and Zeiss Microscopes	

