

[See all 16 Products in Family](#)

# Olympus UPLXAPO60XW 60X Water Immersion Objective

See More by [Olympus](#)



Olympus UPLXAPO 60X Water Immersion Objective

Stock **#74-022** NEW **1 In Stock**

⊖ 1 ⊕ £7,344.<sup>00</sup>

**ADD TO CART**

Volume Pricing	
Qty 1+	£7,344.00 each
Need More?	<a href="#">Request Quote</a>

ⓘ Prices shown are exclusive of VAT/local taxes

## Product Downloads

### General

UPLXAPO60XW **Model Number:**

Compatible Tube Lens Focal Length (mm):  
Focal Length: 180mm

Microscope Objective **Type:**

Infinity Corrected **Style:**

Olympus **Manufacturer:**

Spring-loaded. **Note:**

## Physical & Mechanical Properties

31.50 **Diameter (mm):**

0.44 **Field of View (mm):**

49.51 **Length (mm):**

44.61 **Length excluding Threads (mm):**

31.5 **Maximum Diameter (mm):**

7.20 **Exit Pupil Diameter (mm):**

## Optical Properties

0.17 **Compatible Cover Glass Thickness (mm):**

3.00 **Focal Length FL (mm):**

60X **Magnification:**

1.20 **Numerical Aperture NA:**

0.28 **Resolving Power (µm):**

0.25 **Depth of Field (µm):**

0.28 **Working Distance (mm):**

400 - 1000 **Wavelength Range (nm):**

26.5 **Field Number (mm):**

45.00 **Parfocal Length (mm):**

Water **Immersion Liquid:**

## Threading & Mounting

RMS / 20.32mm x 36 TPI **Mounting Threads:**

## Regulatory Compliance

[View](#) **Certificate of Conformance:**

[Contains SVHC\(s\)](#) **REACH 241:**

## Product Details


- High Numerical Aperture (NA) up to 1.45
- Chromatic Aberration Correction from 400 - 1000nm
- Uniform Image Flatness over Large FOVs

Olympus X-Line Extended Apochromat Infinity Corrected Objectives (also referred to as the UPLXAPO series) are Olympus' next generation high performance objectives that simultaneously feature high numerical aperture (NA), broadband chromatic aberration correction, and excellent image flatness. The high numerical aperture (NA) allows for excellent light collection efficiency to better capture faint signals, enabling a shorter exposure time and higher contrast. The broad chromatic aberration correction compensates for the wavelength-dependent spatial shift of colors to achieve excellent color reproduction. These objectives also achieve high image flatness from the center to the edges of images over large fields of view (FOV), making them ideal for image stitching and quantitative image analysis. Olympus X-Line Extended Apochromat Infinity Corrected Objectives are used in demanding life science and clinical applications including pathology, fluorescence microscopy, multicolor or widefield fluorescence imaging, confocal or super resolution microscopy, and whole slide imaging.

## Technical Information

	Stock #	A	B	C	D	E	F	G	WD
Olympus UPLXAPO 4X Objective	<a href="#">#14-905</a>	31.89	-	-	15.4	15.4	28	4.7	13.0
Olympus UPLXAPO 10X Objective	<a href="#">#14-906</a>	41.78	3.7	4.28	9.6	19.4	28	4.6	3.1
Olympus UPLXAPO 20X Objective	<a href="#">#14-907</a>	44.28	3.39	2.88	10.2	18.6	28	4.8	0.61

Olympus UPLXAPO 40X Objective	#14-908	44.71	4.21	2.11	7.5	20.3	31.5	4.9	0.18
Olympus UPLXAPO 60X Oil Immersion Objective	#14-909	44.74	2.1	2.84	10	19.3	31.5	4.9	0.15
Olympus UPLXAPO 100X Oil Immersion Objective	#14-910	44.76	2.1	3.06	7.5	19.1	30.5	4.8	0.135

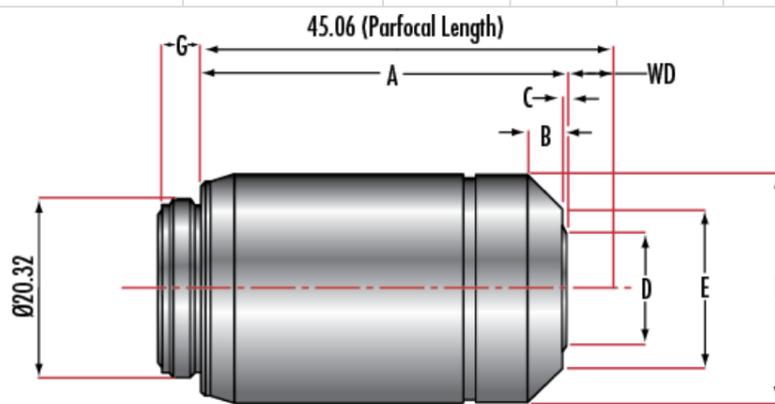

;